

OST Medical

An Enteral Technology Company

The Sentinel Enteral Feeding Pump from OST Medical

- Dose Check keeps a record of the accumulated amount of formula delivered over several feedings.
- Clear, Large Display Windows can be read from a distance. Displays are adjustable for brightness for day or night operation.
- The Clear Button clears all displays, except total accumulated dose, with one, four-second touch.
- Audible and Visible Alarms are readily differentiated. User is alerted when battery is low, when the formula container is empty, or the line is occluded, when the dose is complete, and in the event of formula excess flow.
- The Smooth, Plastic Case is easily cleaned with mild soap and water.
- Built-in Pole Clamp allows patient mobility.
- Pump is easily loaded with one hand.
- Memory is not erased in a power failure or when unit is turned off.
- Dose Limit can be set in 1 ml increments up to 2000 mls, then 5 ml increments up to 9999 mls.
- Programming Directions are located directly on the pump.
- Special Cap Holder for administration set on back of pump, keeps caps convenient and helps prevent contamination.

Available Feeding Tubes

Part No.	Description	Packaging
E-1001	Spike Set	30 Each
E-1200	1200 ml Bag Set with attached delivery set	30 each
E-0500	500 ml Bag Set with attached delivery set	30 each

SENTINELTM
feeding pump

PATENT
PENDING



Simplicity, Convenience & Economy

OST Medical

OST Medical, Inc.
11 Knight Street, Building F-23
Warwick, RI 02886
Tel.: 401-737-3774
FAX: 401-737-4519
E-mail: sales@ostmedical.com
www.ostmedical.com

ISO 13485

OM-0008 Rev-A


The **SENTINEL**TM
Enteral Feeding Pump
& Delivery Sets

SENTINEL™
feeding pump

The **Sentinel™** is a time-saving, easy-to-use enteral feeding pump that can be used in the administration of all enteral feeding formulas. The pump is easy to program, and many of its functions can be activated with a single touch. Reliable and accurate, The **Sentinel™** is designed to give the user easy access to its many useful features.

Features & Benefits

- Advanced Software - Pump odometer records and analyzes data for usage history
- Battery Life - Internal monitor for battery life estimation
- Large LED Display - sizeable area for easy viewing
- Alarm Identification - differences in tone sequence make alarms easy to identify

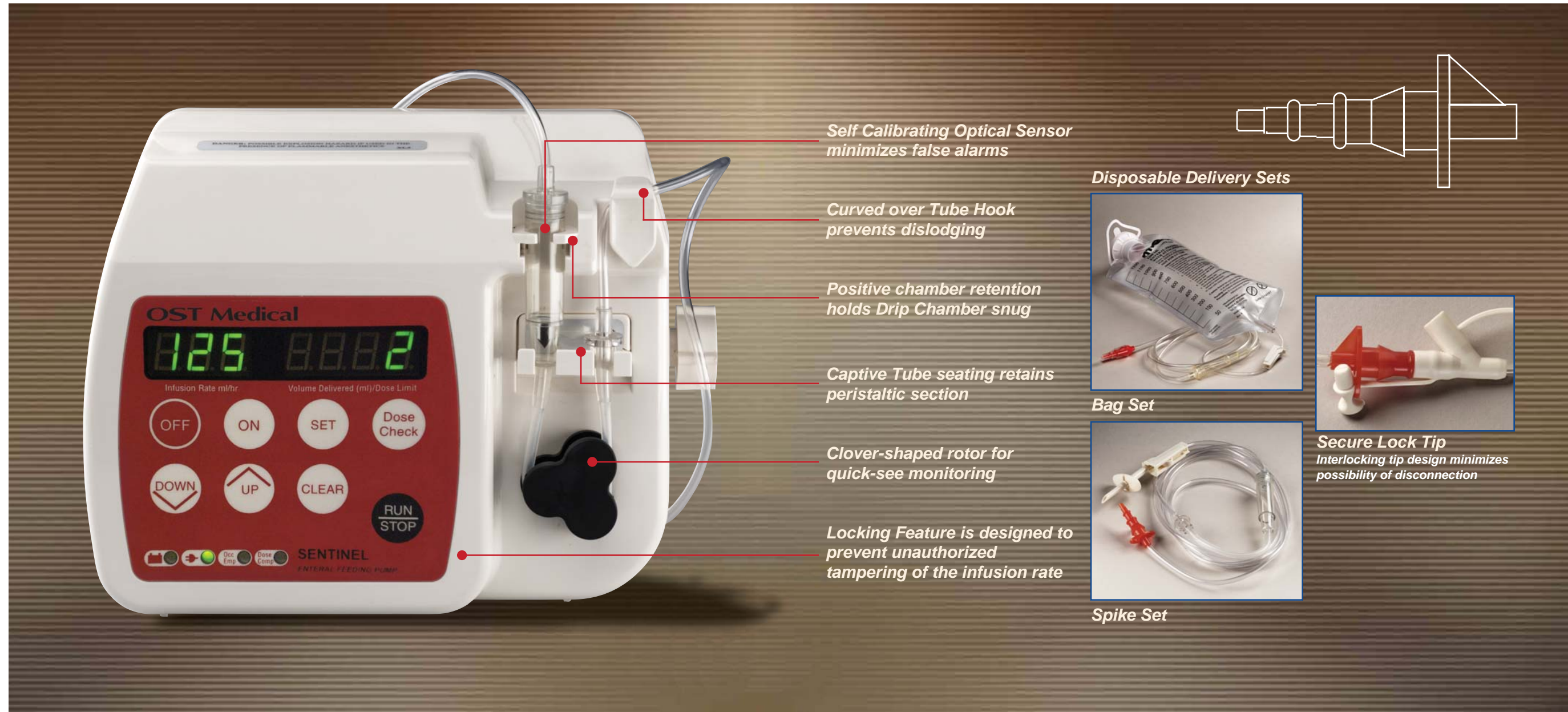

- Easy-to-Use
- Easy-to-Program
- Reliable & Accurate
- Competitive Price

3 Step easy-to-use operation


1. Press ON.
2. Press UP/DOWN to set flow rate.
3. Press RUN/STOP to start/stop feeding.

OST Medical

Impact resistant housing
withstands hard impacts

Integral pole clamp
easy grip design

Cap holder
for convenience and helps prevent contamination

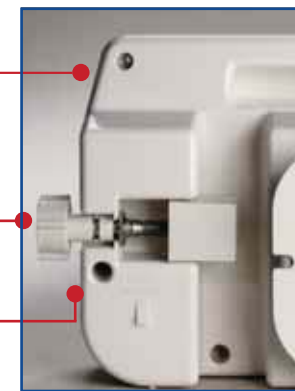

Rates	5-295 ml/hr in 1 ml increments
Dose Limit	1-2000 ml(s) in 1 ml increments/2000-9999 mls in 5 ml increments
Accuracy	+/- 10%
Battery	10 hours at 125 ml/hr
Recharge	Continuous while unit is plugged into AC outlet
Memory	Retained without power required
Size	7.5" x 9.1" x 5.2"
Weight	5.2 lbs.